

Floor Scale

Model SRV95X Series

The ultra low profile SRV95X series are sleek innovative scales, that provide years of reliable service and accurate weight data.

The Scale features a wired, wall mounted large 4.3" color touch screen display that is easy to read.

The low-rise scale surface has a non-slip polyurea coating infused with an aggregate to provide a slip-free surface for animals thus reducing the fear that animals often experience with a scale.

An optional second readout provides ease of viewing for busy staff.



CAPACITY

400 lb or 182 kg

PLATFORM SIZES

- SRV955 24" x 36"
- SRV957 24" x 48"
- SRV959 24" x 60"

AUTO-WEIGH

An auto-weigh feature allows the scale to wake, and display a weight when an animal steps on.

AUTO-HOLD

Scale will automatically lock in on weight above 1 kilogram.

MEMORY RECALL

Recall the last stored weight.

WIRELESS CAPABILITY

Scale comes with peer to peer wi-fi to communicate with an optional auxiliary display (FR95X-AUX). Note: facility does not need wi-fi for this feature to work.

QR WEIGHT DISPLAY

Simply touching the displayed weight will switch the weight to a QR code for easy scanning into your medical record system.

Floor Scale

Model SRV95X Series

SR Scales[®]
by SR Instruments, Inc.

Quality. Reliability. Accuracy.

SYSTEM SPECIFICATIONS


MAX WEIGHT	400 lb / 182 kg
DISPLAY RESOLUTION	0.1 lb / 0.1 kg
BUTTON FUNCTIONS	Zero/Weigh, Hold/Recall
DISPLAY UNITS	Pounds / Kilograms
POWER	120 VAC 50-60 Hz
AUTO HOLD	When enabled scale will automatically lock in a weight when stable, and store it. To recall this weight press the "Recall" button.
ALWAYS ON	Scale is always on in the background and will automatically wake up and lock a stable weight when weight is detected.
WARRANTY	Two (2) Years

ACCESSORIES AND OPTIONS


DESCRIPTION	PART NUMBER
SRV95X WIRELESS AUX DISPLAY	FR95X-AUX



DIMENSIONS in [cm]

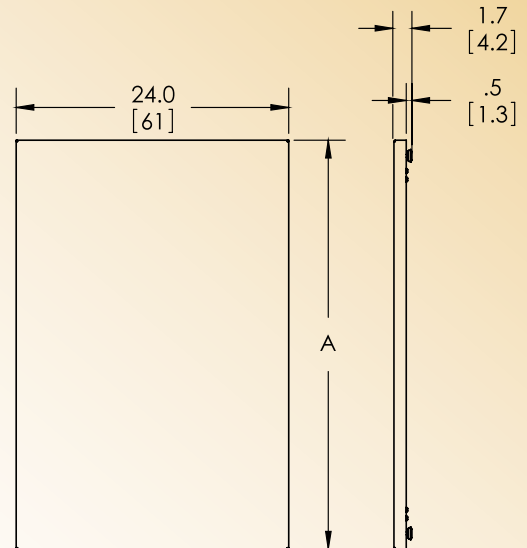

PACKAGING

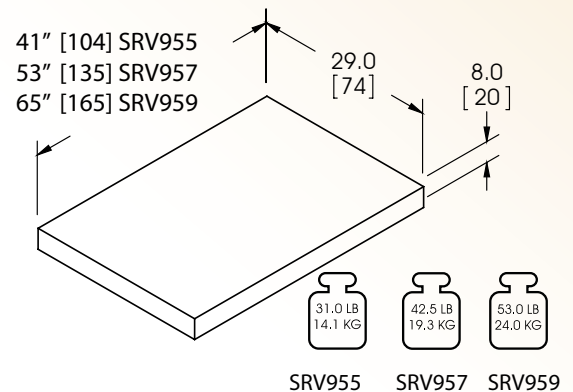

Specifications subject to change without notice.